


RECOMBINANT EQUINE IL-2

KINGFISHER

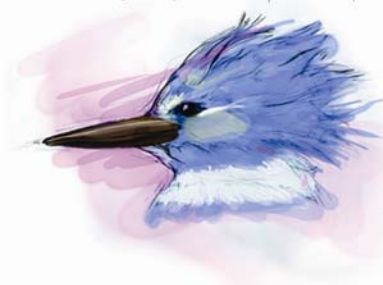

BIOTECH, INC.

Catalog Number: RP0006E-005

Quantity: 5 µg

Source: Yeast

Formulation/Reconstitution: Lyophilized. Reconstitute with sterile phosphate-buffered saline containing at least 0.1% carrier protein.

Stability and Storage: Stable for up to twelve months from date of receipt at -20°C. Stable for at least 3 months when stored in working aliquots with a carrier protein at -20°C. Avoid repeated freeze/thaw cycles.

Molecular Weight: 14.9 kDa (calculated)

Purity: >95% as visualized by SDS-PAGE analysis.

Purification: Ion-exchange chromatography

Entrez Gene ID 100034204

Amino Acid Sequence: APTSSSKRET QQQLKQLQMD LKLLLEGVNN NKNPKLSKML
TFKINMPKKA TELKHLQCLE EELKPLEEML KNFLSKDIKE
LMSNINVTVL GLKGSETRFT CEYDDETGFI VEFLNKWITF
CQSIFSTMT (129)

Bioactivity: The biological activity of equine IL-2 was measured in a cell proliferation assay using carboxyfluorescein diacetate succinimidyl ester-loaded equine PBMC. The ED50 for this effect is typically 0.5 - 1.0 ng/mL.

Country of Origin: USA

Support: Generation of this material is based upon work supported by the Cooperative State Research, Extension, Education Services (CSREES), USDA under project #0206006 US-Veterinary Immune Reagent Network (www.vetimm.org).

Warranty: Products are warranted by Kingfisher Biotech, Inc. to meet stated product specifications and to conform to label descriptions when used, handled and stored according to instructions. Unless otherwise stated, this warranty is limited to one year from date of sale or expiry date, whichever comes first. Kingfisher Biotech's sole liability for the product is limited to replacement of the product or refund of the purchase price. Kingfisher Biotech products are supplied for research applications. They are not intended for medicinal, diagnostic or therapeutic use. The products may not be resold, modified for resale or used to manufacture commercial products without prior written approval from Kingfisher Biotech, Inc.

REAGENTS FOR VETERINARY RESEARCH

